




# Gastric Sump Tube

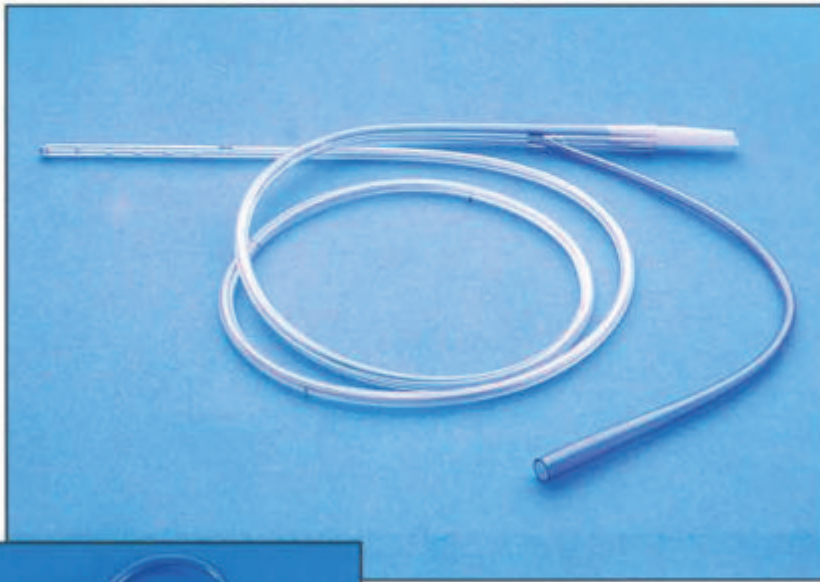

Blue side arm caps to 5-in-1 connector when not in use, reducing reflux accidents

- Double lumen PVC tube - a large lumen for suction and a small lumen for irrigation and venting
- Blue side arm extends into the smaller irrigation/vent lumen
- Blue side arm allows simultaneous irrigation of stomach without discontinuing suction and doubles as a vent lumen when irrigation is not taking place
- Unique double lumen design and vent feature enhances patient safety by eliminating risk of injury to stomach lining
- Funnel end with removable 5-in-1 connector easily connects to suction units
- Blue side arm caps to 5-in-1 connector when not in use, reducing reflux accidents
- Radiopaque stripe facilitates proper placement under x-ray
- Adult and pediatric sizes
- Sterile

| Part Number | French Size | O.D. (mm) | Length     | Qty per cs |
|-------------|-------------|-----------|------------|------------|
| 1217-10     | 10          | 3.3       | 36"(91cm)  | 50         |
| 1217-12     | 12          | 4.0       | 48"(122cm) | 50         |
| 1217-14     | 14          | 4.7       | 48"(122cm) | 50         |
| 1217-16     | 16          | 5.3       | 48"(122cm) | 50         |
| 1217-18     | 18          | 6.0       | 48"(122cm) | 50         |